MELAstore Tray Ophthalmologie / MELAstore Tray Ophthalmology

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Instructions for use

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
MELAstore Tray Ophthalmologie / MELAstore Tray Ophthalmology

Instructions for use

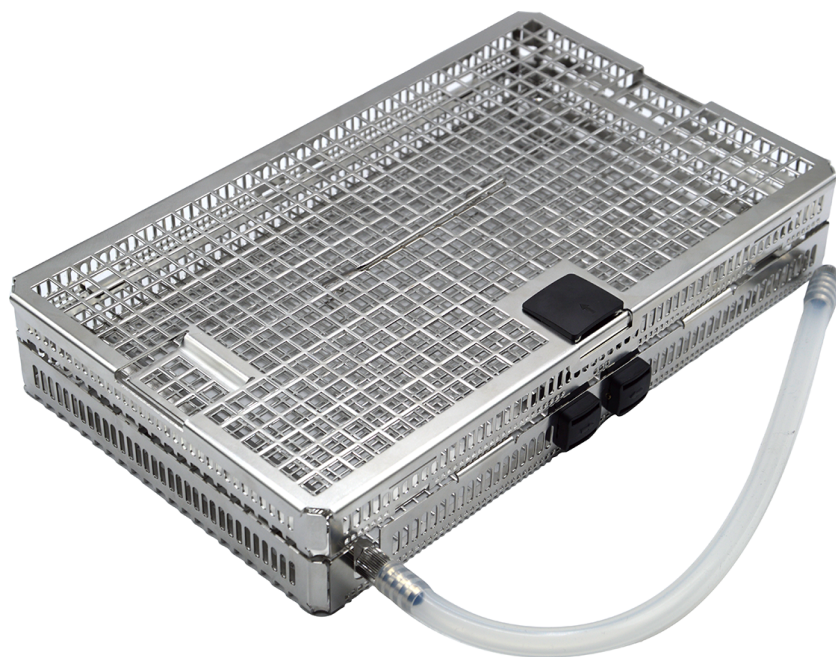Email: info@melag.com
Web: www.melag.com

Original instructions

Responsible for content: MELAG Medizintechnik GmbH & Co. KG
We reserve the right to technical alterations

Intended use

The MELAstore Tray Ophthalmology is used for reprocessing thermostable solid instruments and hollow-body instruments for cataract surgery in the MELAtherm 10 and MELAtherm 10 Evolution. After reprocessing in the MELAtherm, the MELAstore Tray Ophthalmology, with instruments, can be sterilized, stored and transported using a MELAstore Box 200.

Up to six MELAstore Tray Ophthalmology can be reprocessed together in the MELAtherm 10 and MELAtherm 10 Evolution.

Scope of delivery

Please check the scope of delivery before using the product.

Description		Description	
2x Luer connection (male)		4x hose connector (4 mm) with external thread	
3x Luer connection (female)		Silicone hose (10/6 mm) ^{*)}	
3x Luer-Lock adapter (male) ^{*)}		Silicone hose (internal) 500 mm	
Hose connector (6 mm) with external thread ^{*)}		Instrument holder for flex baskets (6 pcs.) ^{*)}	
Hose connector (6 mm) with internal thread ^{*)}		4x silicone bar	

^{*)} Standard accessory for MELAtherm 10/MELAtherm 10 Evolution, see user manual Accessories for MELAtherm

Safety instructions

Qualified personnel

- Instrument reprocessing using a washer-disinfector may only be carried out by competent personnel.

Accessories

- Use original accessories from MELAG. MELAG cannot provide a guarantee for non-MELAG accessories, even if they are in possession of validation.
- When using additional accessories from other manufacturers to secure instruments, especially hollow-body instruments, follow the instructions in the operating manual provided by the manufacturer of the accessories.
- If using accessories from other manufacturers, please ensure that they are heat-resistant up to 95 °C and resistant to the process agents used.

Personal protective equipment

- For your own safety and to avoid injuries, always wear personal protective equipment (PPE) or other suitable hand protection when loading and unloading instrument and washing baskets.
- Take care when arranging sharp and pointed instruments and do so in a manner that prevents all danger of injury. Loading should best be performed from back to front and unloading from front to back. Wear suitable protective gloves.
- When spraying the instruments with the MELAjet spray pistol, make sure that the adapter is connected correctly and that the instruments to be sprayed are used properly, see the MELAjet user manual. Wear protective gloves and goggles.

Reprocessing

- Only reprocess instruments that are suitable for automatic reprocessing.
- In Germany, the DGSV (German Society for Sterile Supply) recommends pre-cleaning contaminated instruments after an operation in accordance with AK-Q-70 and AK-Q-71.
- If required, order the reprocessing instructions from the manufacturer in accordance with EN ISO 17664.
- Reprocessing in the washer-disinfector may only be carried out using an injector rail with a plastic central filter.
- Observe the instructions provided in the washer-disinfector user manual regarding ophthalmological instruments.

Loading instructions

- Prepare instruments with joints that cannot be dismantled or instruments that can be closed (e.g. tweezers) in an open state.
- Only reprocess instruments in the washer-disinfector that have been approved by the manufacturer. As a rule, these instruments are marked directly or in the reprocessing instructions with the following symbol: 
- Comply with the reprocessing instructions provided by the instrument manufacturer, especially those pertaining to compatibility with process agents.
- Do not reprocess any disposable instruments. Disposable instruments are usually marked with the following symbol: 


Application

Note the following:

- Install the hoses free of kinks and sacks.
- Keep the hoses as short as possible. Water can accumulate in long hoses.
- Do not close hoses.
- Check the hoses, connections and instruments before and after reprocessing to ensure that they are tight. If a hose, connection or instrument has come loose, reprocess the instrument again.
- Push hoses onto the hose connections as far as they will go so that the hoses do not come loose during reprocessing.
- Instruments may only be reprocessed fixed in the trays. Use either instrument holders or silicone bars for this purpose.
- Only insert as many instrument holders and silicone bars into a tray as are required for one reprocessing process.

The MELAstore Tray Ophthalmologie consists of a lid, an upper basket and a lower basket including an injector rail.

For reprocessing in the washer-disinfector, the upper and lower baskets are loaded separately and placed in the washing chamber.

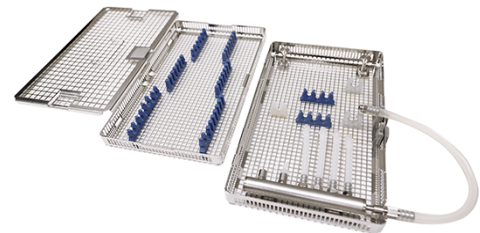
Separating the upper and lower baskets

The upper and lower baskets of the MELAstore Tray Ophthalmologie are assembled for delivery. For reprocessing in the MELAtherm, the upper and lower baskets must be separated from one another.

The upper basket (pos. a) is suspended in the lower basket (pos. b).

PLEASE NOTE: The illustration shows the upper and lower baskets without a lid.


1. Press the lock on the lower basket (pos. c) together and separate the upper basket (pos. a) from the lock.
2. Unhook the upper basket (pos. a) from the lower basket.

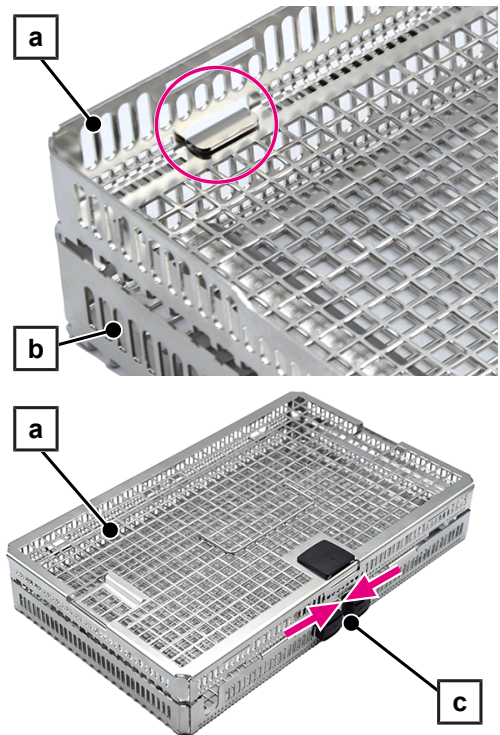Equipping and loading the upper basket


The following must be fulfilled or present:

- ✓ *The upper and lower baskets are separated.*
- 1. Slide the lid lock (pos. a) in the direction of the arrow and remove the lid from the lock.
- 2. Unhook the lid from the upper basket (see circular marking).

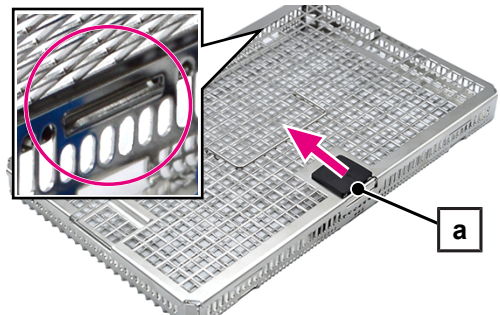
3. Insert the silicone bars and/or instrument holders. If necessary, cut off individual silicone bars and insert them separately.
4. Press the feet (pos. b) firmly into the base of the upper basket so that they protrude from the underside.

PLEASE NOTE: The silicone bars/instrument holders must be flush with the base.

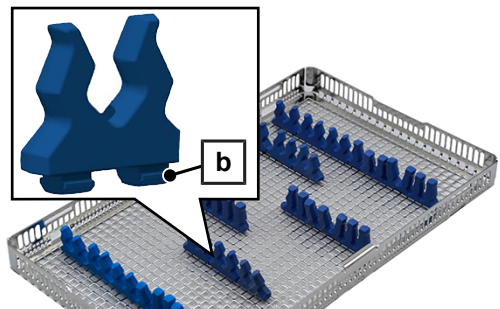5. Load the upper basket with solid instruments correctly for rinsing (e.g. hinged instruments opened, no loading over each other).
 6. Hook the lid (pos. c) into the upper basket (see circular marking) on both sides.
 7. Close the upper basket with the lid (pos. c). Press the lock (pos. a) of the lid onto the upper basket.
- A clicking sound can be clearly heard when it is correctly closed.


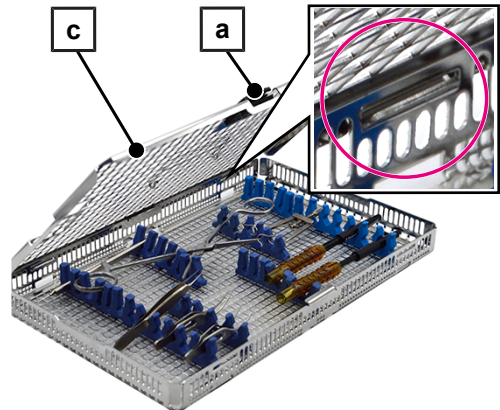
Equipping and loading the lower basket

Equipping the lower basket

The following must be fulfilled or present:

- ✓ *The upper and lower baskets are separated.*
1. Screw the required number of Luer-Lock adapters (male, pos. c) into the injector rail.
 2. Screw the required number of hose connectors (4 mm, pos. b) into the injector rail.
 3. Cut the silicone hose from the scope of delivery (500 mm) to the required lengths.
 4. Push one cut silicone hose (pos. a) onto a hose connector (pos. b) on the injector rail as far as it will go.
 5. Equip the free end of each silicone hose with the required plastic Luer connection (pos. d).
 6. Insert the silicone bars and/or instrument holders. If necessary, cut off individual silicone bars and insert them separately.
 7. Press the feet (pos. e) firmly into the base of the lower basket so that they protrude from the underside.
- PLEASE NOTE:** The silicone bars/instrument holders must be flush with the base.

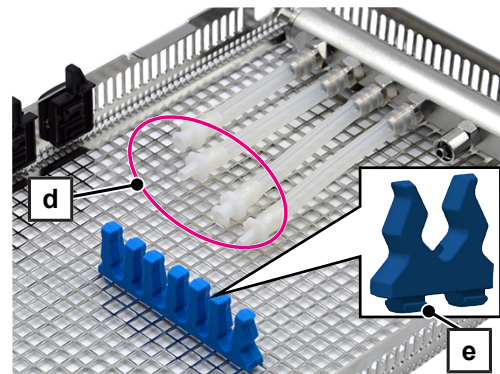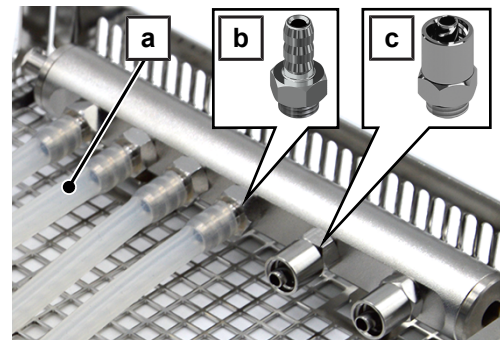
Loading the lower basket


If possible, reprocess the Phaco needle of the Phaco handpiece separately.

1. Unscrew the Phaco needle (pos. f) from the Phaco handpiece.
2. Insert the Phaco needle (pos. f) into the Luer connection (female, pos. g).
3. Connect the Luer connection (female, pos. g) to the Luer-Lock adapter (male, pos. c) on the injector rail.

PLEASE NOTE: Due to different geometries of the Phaco needle, it may not be possible to insert it fully into the Luer connection (female). In this case, the Luer connection cannot be connected to the Luer-Lock adapter (male). Leave the Phaco needle screwed onto the Phaco handpiece for reprocessing.

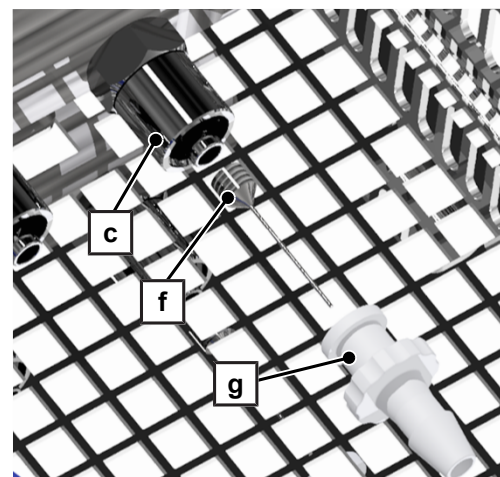4. Connect the appropriate hollow-body instrument with a half turn with the Luer-Lock adapter (male, pos. i) onto the injector rail (pos. h).
5. Connect the appropriate hollow-body instrument with the Luer connection (pos. j) to the silicone hose.
6. Fix the instruments in the silicone bars/instrument holders.
7. Install the cable of the Phaco handpiece on the outer edge of the lower basket around the injector rail (pos. h) and bars (pos. k).
8. Seal the unused connections on the injector rail with a screw plug (art. no. ME80140).
9. If available, close the Luer-Lock adapter (male) with the closure (female) for Luer/Luer-Lock (art. no. ME80180).

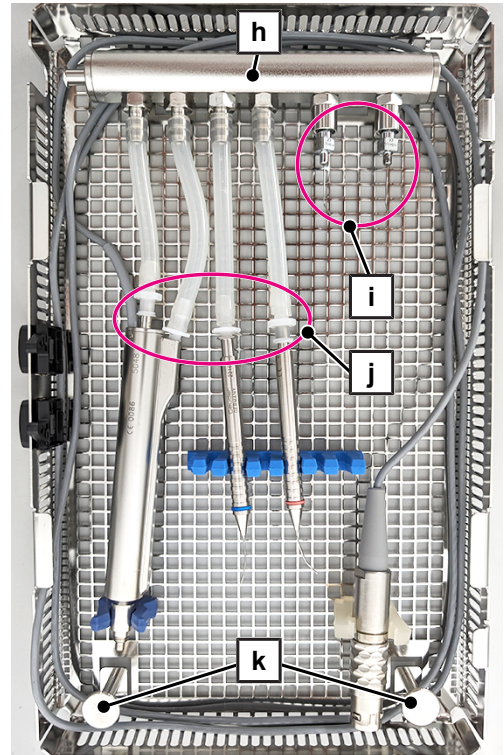Inserting the upper and lower baskets into the washing chamber

Loads are inserted into the washing chamber using the universal holder Flex 3, see user manual Accessories for MELAtherm. Up to six upper baskets or lower baskets can be placed per holder. A maximum of six lower baskets per batch can be used for reprocessing.

WARNING


Danger of contamination

Place the closed upper basket and the open lower basket of the MELAstore Tray Ophthalmologie separately in the MELAtherm basis basket.


The following must be fulfilled or present:

- ✓ *The injector rail of the basis basket is equipped with a plastic central filter.*


1. Screw the hose connector with external thread (6 mm) into the injector rail of the basis basket by hand and slide the silicone hose (10/6 mm) onto the connector as far as it will go.

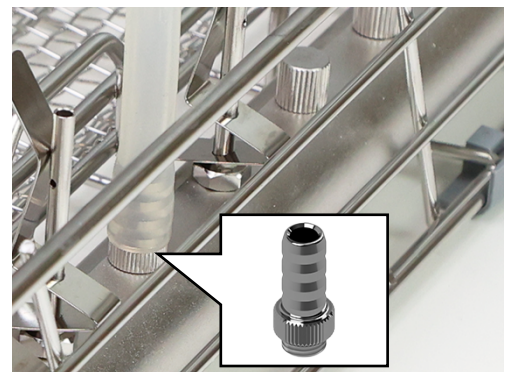2. Screw the hose connector with internal thread (6 mm, pos. b) onto the Luer-Lock adapter (male, pos. a) by hand.

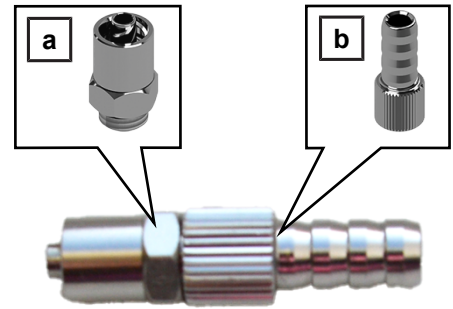3. Slide the free end of the silicone hose (10/6 mm) onto the hose connector (pos. b) as far as it will go.

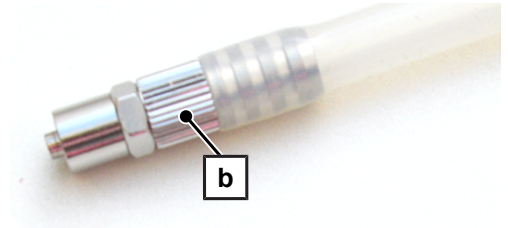
4. Place the universal holder Flex 3 next to the injector rail (pos. c) of the basis basket.
5. Place the **open lower basket** in the holder directly next to the injector rail (pos. c) of the basis basket. For easier loading and unloading, the lower basket should be inclined towards the device door.

- The injector rail of the lower basket is at the top.
- The connection side of the injector rail of the lower basket is fixed in the direction of the injector rail (pos. d) of the basis basket.

6. Connect the Luer-Lock adapter (male, pos. a) to the injector rail (pos. d) of the lower basket with a half turn.

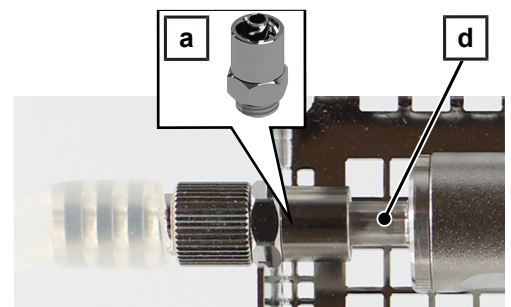7. Place the **closed upper basket** in the mount. For easier loading and unloading, the upper basket should be inclined towards the device door.



Following reprocessing in MELAtherm

1. Check adapters, hoses, connections and instruments to ensure that they are firmly seated. If an adapter, hose, connection or instrument has come loose, the respective instrument must be reprocessed again.
2. Check that all instruments have been successfully cleaned.
3. Check hollow bodies for free passage. Observe the department-specific instructions in the user manual of the washer-disinfector.
4. Dry the instruments, particularly complex hollow bodies such as Phaco handpieces using medical compressed air to remove residual moisture.


Use with MELAstore Box 200

The upper and lower baskets of the MELAstore Tray Ophthalmology must be assembled for use in the MELAstore Box 200.

1. Remove the Luer-Lock adapter (male) incl. connector and silicone hose from the injector rail of the lower basket with a half turn.
2. Hook the closed upper basket (pos. a) into the lower basket (pos. b) on both sides.

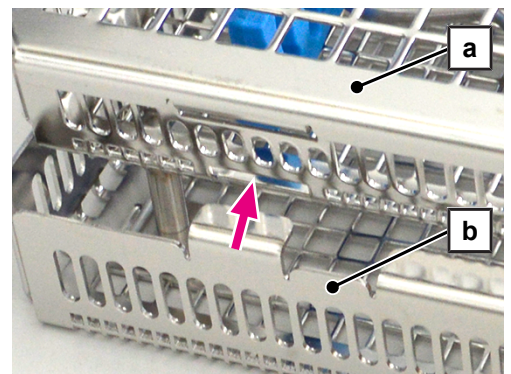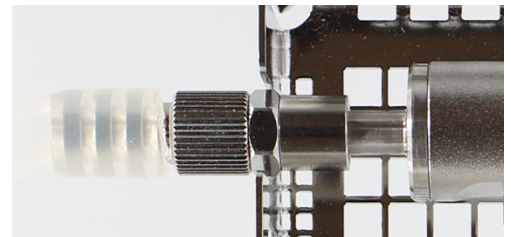- Press the upper basket (pos. a) down until it engages in the lock (pos. c) of the lower basket (pos. b).

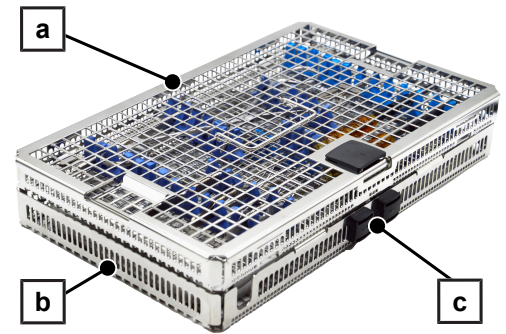- Insert the MELAstore Tray Ophthalmologie into the MELAstore Box 200.
- Comply with the application usage advice for the MELAstore Box 200 according to the separate document.

Routine checks

Regularly check the plastic parts for wear and replace them if necessary.

Regularly check all accessories for damage. If necessary, dispose of damaged accessories properly.

Storage

Store all accessories in a dry atmosphere protected from corrosive objects or media.

Store the spare parts in a dry and dust-protected place, e.g. in the packaging.

Disposal

Dispose of accessories if they show signs of wear (e.g. abrasion, cracks or corrosion) or damage (e.g. deformation) to a great degree.

Reprocess the accessories to be disposed of in isolation in the washer-disinfector before disposing of them.

Accessories and spare parts

You can obtain the specified articles together with an overview of further accessories from your stockist. Information regarding the instrument reprocessing accessories can be found in the current MELAG price list.

Accessories	Art. no.
Luer-Lock adapter (male)	ME74130
Universal holder Flex 3	ME80136
Screw plug for injector rail and distributors	ME80140
Hose connector (6 mm) with external thread	ME80150
Hose connector (6 mm) with internal thread	ME80160
Closure (female) for Luer/Luer-Lock	ME80180
Silicone hose (10/6 mm), 2 m	ME80190
Instrument holder for Flex basket (60 pcs.)	ME80395
Basis basket with injector rail incl. 2 plastic central filters	ME80440

Only original MELAG spare parts may be used.

Spare parts	Art. no.
Lid handle	ME12423
Luer connection (male)	ME12431
Luer connection (female)	ME12432
Silicone hose (internal) 500 mm	ME12435
Locking set	ME22999

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